



Wood Frame Science Table Instructions

939300-INST
Rev. A

- 1.** Place table top on the ground with the frame facing upwards. *If installing leg protector boots, complete INST-BOOT before completing these instructions.*

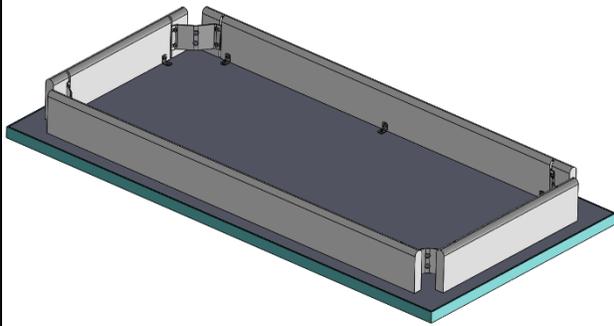


- 2.** Unscrew the nut and washer from each leg bolt. Align each leg in position on the table top, inserting the leg bolts into the corner brackets.

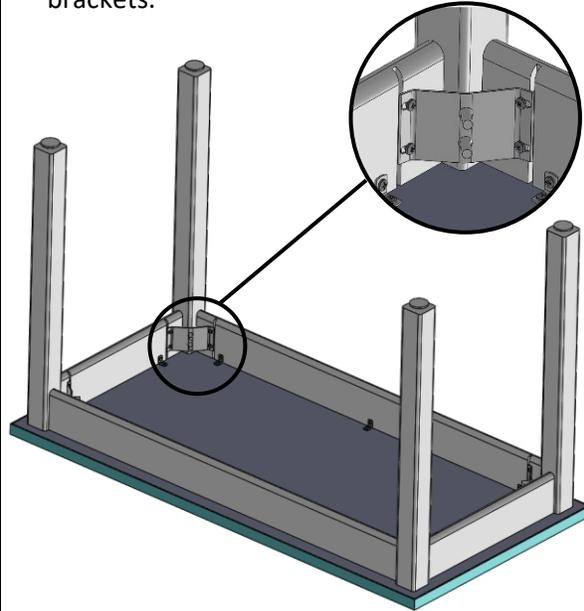

- 3.** Attach the legs to the frame with (2) nuts and (2) washers on the (2) leg bolts in each corner bracket. Tighten to a torque of 50-60 in-lbs.



- 4.** Flip table to in-use position. *Your product may vary from what is pictured.*

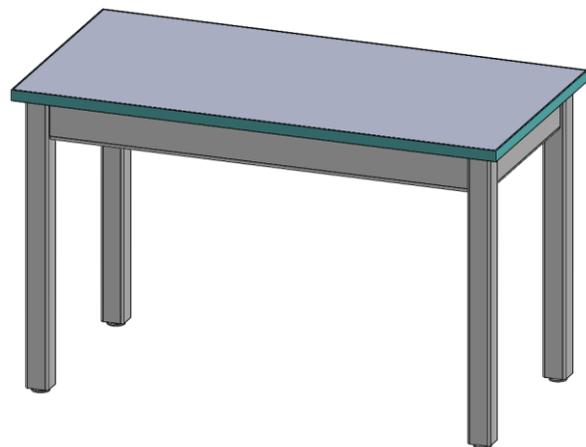

For questions please call us at (616) 396-1142